FERROVAC GMBH

ULTRA HIGH VACUUM TECHNOLOGY

SALES@FERROVAC.COM WWW.FERROVAC.COM


рноме +41 44 273 16 38

WOBBLESTICK MANIPULATORS

Ferrovac wobblestick manipulators define the benchmark for delicate and precise manipulations in ultra high vacuum systems.

Due to the principle of magnetic force coupling, the motion of the shaft is not affected by the difference between atmospheric pressure and ultra high vacuum. The operator can guide the shaft of the wobblestick nearly force-free in a smooth and controlled way.

Magnetically driven wobblestick manipulators always allow for continuous rotation of the shaft.

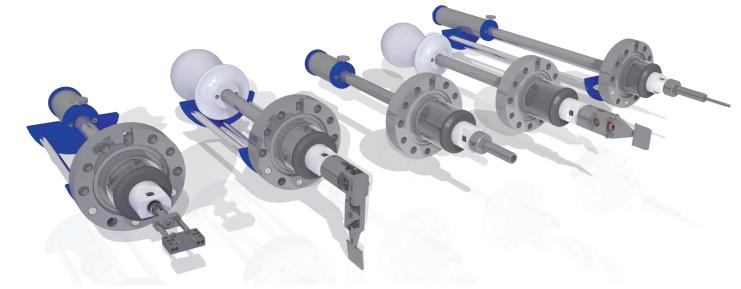

WMG40 Dual Shaft Wobblestick

If you would like to grip a sample holder with a pincer, a dual shaft wobblestick is needed. Slight pressing of the two halfs of the knob against each other, will open the pincer.

A wide selection of pincers is available for any known type of sample holder.

Ask for tailored jaws if your application is unique.

P Starrison

FERROVAC GMBH

ULTRA HIGH VACUUM TECHNOLOGY

THURGAUERSTRASSE 72 CH-8050 ZÜRICH

SALES@FERROVAC.COM WWW.FERROVAC.COM

PHONE +41 44 273 16 38

